


Introducing The FORS™-15 Off-Loading Insole¹

The Durable and Adaptable Foot Off-Loading Solution



For Surgical Shoes
(Shoe not included)

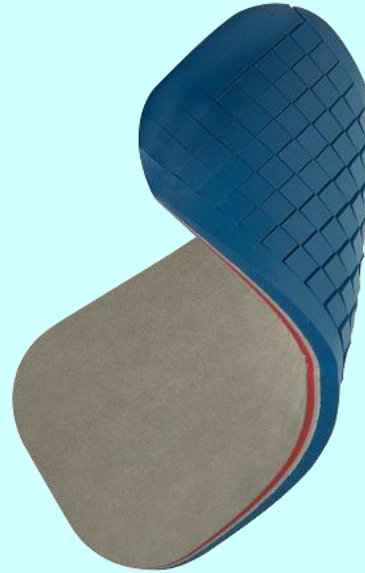

For Therapeutic Boots
(Boot not included)



Consistent Off-Loading & Cushioning

- Durable memory resistant Poron® for consistent cushioning, crush resistance, and long life.
- Highly stable off-loading zones; pre-cut plugs are torn out to create off-loading zones, but remain solidly in place when not disturbed. If you have ever experienced products with unstable and shifting pegs that produce inconsistent off-loading, you will welcome the durability and predictability of the FORS™-15 system!

Best in Class Materials

- Total of 15mm thick Poron® padding, with 5 mm always present under the foot, even after plugs are removed. Unique fabric mid-layer prevents off-loaded areas from collapsing. Tough and abrasion resistant Alcantara® top covers minimize slip and shear, absorb perspiration, are non-allergenic and bacteria resistant.
- Multilayered Poron® provides softness and support where needed. Softer "slow rebound" Poron® adjacent to the foot, and denser Poron® below for support.

Uniquely adaptable:

- The FORS™-15 insole can be trimmed to any desired shape or size.

Easy Care:

- To completely clean and disinfect, soak the FORS™-15 insole in 10% chlorine bleach solution for 10 minutes and lightly scrub as needed, then thoroughly rinse and air dry.




www.bonapeda.com
(973) 415-9090



BONAPEDA Enterprises LLC
Exclusive U.S. Distributor, Orthotics by Saluber

The FORS™-15 Off-Loading Insole¹ - Construction

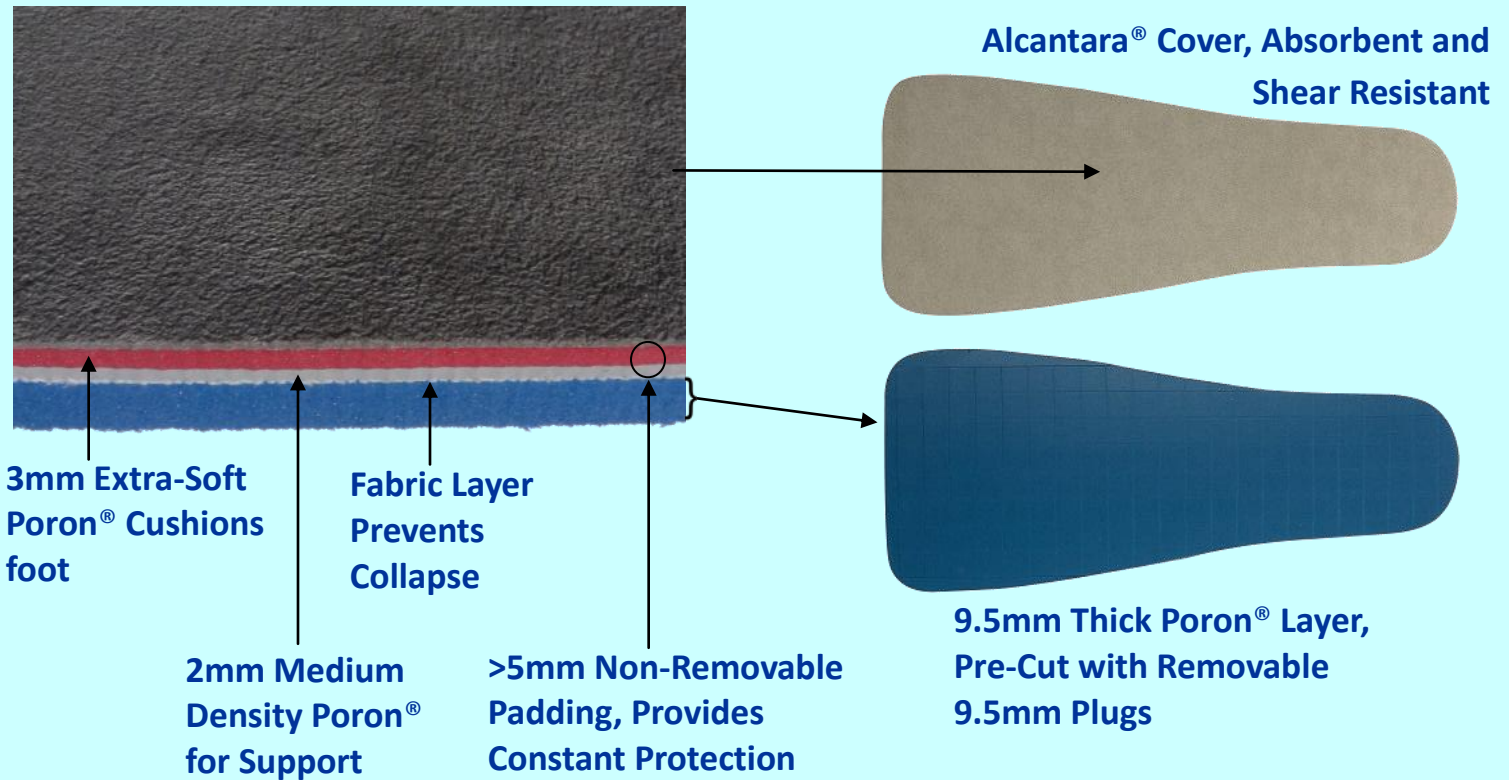

PRODUCT #	DESCRIPTION	SIZES ²	SHAPE
M494-15PZ	FORSTM-15 Insole: 15mm Multi-Layered Poron® Insole with Removable Plugs, Alcantara® Top Cover (Sized to fit medical/surgical shoes)	Women's Medium/Large ³ Men's Small Men's Medium Men's Large	
M494-15PZ	FORSTM-15 Insole: 15mm Multi-Layered Poron® Insole with Removable Plugs, Alcantara® Top Cover	XXL (BOOT SIZE)	

1. Patent Pending
2. FORSTM-15 sizes are designed to correspond to industry standard surgical shoes and boots constructed with removable foam inserts. It is recommended to remove and replace the original foam insert with the FORSTM-15 Insole. It is not recommended to place the FORSTM-15 Insole on top of the original insole, as it may raise the heel out of the shoe or boot and create instability.
3. The Women's Medium size is created by trimming the heel of the Women's Large size. A trimming guide is printed on the insole.



www.bonapeda.com
(973) 415-9090



BONAPEDA Enterprises LLC
Exclusive U.S. Distributor, Orthotics by Saluber